


# HepaStop



19.695.00

HepaStop Cannula Holder  
by Dr. Freude

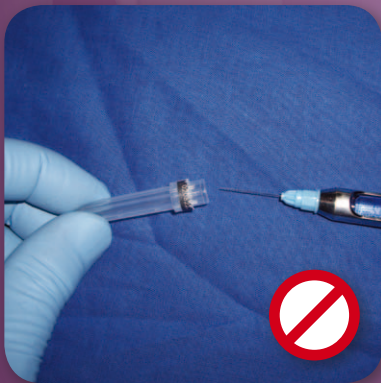

## HepaStop

BY DR. FREUDE / HELMUT ZEPF

**It's quite simple to increase safety in your practice!**

**HepaStop** from **HELMUT ZEPF** meets all requirements for a safe handling with Cartridge Syringes while dental interventions.

Applicable to all basic syringes and tubes.

**HepaStop** is light, handy, space-saving, disinfectable and sterilizable.

Safe working also with fractioned injections with repeated Re-capping during the whole duration of the surgery.

Conform to TRBA 250.

# NEWS

03/11



# ZEPF Endo-Control

## Mouth Mirror Handle Endo-Control

by Dr. Carsten Franke

With the **HELMUT ZEPF** Endo-Control mouth mirror handle it is possible to simply determine the necessary working length due to the laser-marked measuring scale.



The measuring precision can be adjusted to 0,5 mm.

The sandblasted surface reduces reflections on the instrument surface.

### 26.180.07

Universal handle, **ZEPF**-Design, single ended, M 2.5, with Endo calibration, by Dr. Carsten Franke



Application **ZEPF Endo-Control**

## ZEPF Lipcare Mouth Mirror

### 24.062.24

Lipcare Mouth Mirror, diam. 24 mm, Titanium Front Surface\*, sterilizable by Dr. Preusse



\*Front Surface - These mirrors have their reflective coatings on their front surfaces, yielding undistorted imaging, with no disturbing double images.



**HELMUT ZEPF**  
MEDIZINTECHNIK GMBH