

ZEPF Bone Scraper II

In oral surgery, bone replacement material and bone regeneration material is often combined with autogenous bone chips in augmentation.

The Bone Scraper is an ideal instrument for easy and fast harvesting of autogenous bone material. It is used to scrape, collect and transplant autogenous bone.

By using the Bone Scraper, instruments like bone filters, trephines, saws and bone mills are no longer necessary. Its blades allow a harvesting of cortical and cancellous bone.

During scraping, the bone is collected in a chamber and at the same time mixed with blood. The bone material can be implanted directly from the previously sterilized Bone Scraper. The cost-efficient disposable blade offers the practitioner an always sharp and economical instrument.

Essential advantages:

- ergonomisch perfekt
- ergonomically perfect
- easy cleaning
- quickly replaceable disposable blade
- simple application
- cost-efficient



Bone Scraper II

47.957.15

angled version

ZEPF-Design, with sleeve and disposable blade, scale 1:1

47.957.10

straight version

Spare parts:

Disposable blade

47.957.30

PU 3 pieces

47.957.50

PU 5 pieces

Sleeve

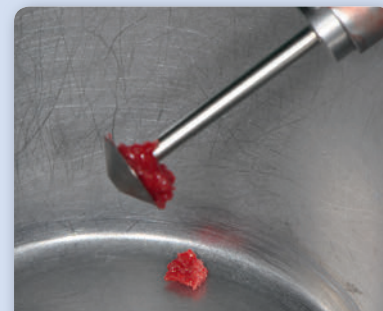
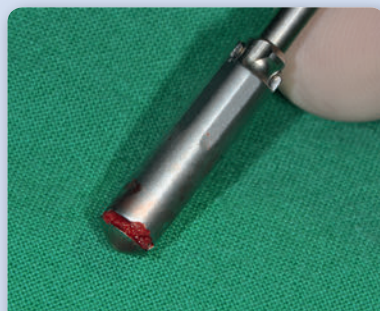
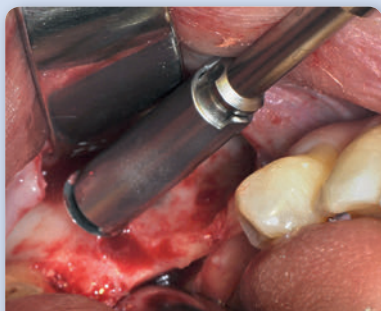
47.957.70

Handle ZEPF-Design

47.957.60

straight version

47.957.65

 angled version,
without illustration


Illustration

Article Description

Order Quantity



47.957.10

1 piece

Bone Scraper II straight, Ø 6 mm with sleeve and one blade, **ZEPF-Design**



47.957.15

1 piece

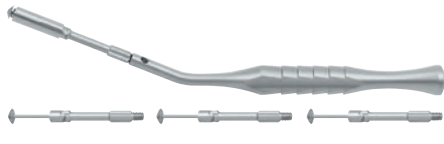
Bone-Scraper II angled, Ø 6 mm with sleeve and one blade, **ZEPF-Design**



47.957.20

1 set

Bone-Scraper II straight, Ø 6 mm, with sleeve and blade Ø 6 mm, including 3 spare blades



47.957.25

1 set

Bone-Scraper II angled, Ø 6 mm, with sleeve and blade Ø 6 mm, including 3 spare blades

Spare parts for Bone Scraper II



47.957.30

1 set

3 Spare Blades, round, Ø 6 mm, plus 1 spare sleeve



47.957.50

1 set

5 Spare Blades, round, Ø 6 mm, plus 1 spare sleeve



47.957.60

1 piece

Handle straight, **ZEPF-Design**



47.957.65

1 piece

Handle angled, **ZEPF-Design**



47.957.70

1 piece

Sleeve for Bone Scraper II

